





FU/363

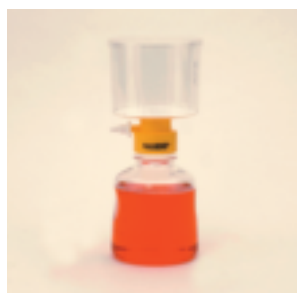

FU/364

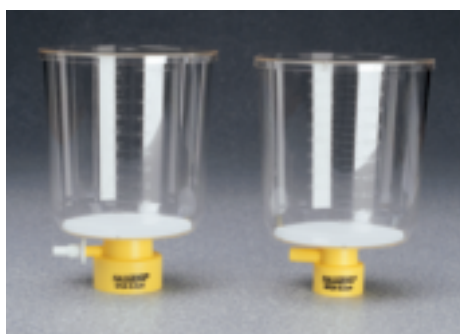

FB/363



FM/362

Filter units, SFCA, Media-Plus

NALGENE

These 90mm diameter (57.2cm² usable filter area) SFCA membrane filter units allow easy vacuum filtration of up to 1 litre of medium with 20% serum. Large membrane surface area and low protein binding provide excellent flow rates and throughput. Upper and receiver 1000ml. In polystyrene housing with polystyrene receiver connected by 1½ turn screw closure. Side-arm features a cellulosic vent plug and quick-disconnect tubing adapter. Colour coded yellow. Supplied sterile in cases of 12 units.

Cat. No.	Code	Pore size, µm
FU/36308	161-0020	0.2
FU/36318	161-0045	0.45

Filter units, SFCA, serum

NALGENE

These 90mm diameter (57.2cm² usable filter area) SFCA membrane filter units are designed for the filtration of large volumes (over 200ml) of whole serum. They provide excellent flow rates and throughput when filtering up to 500ml of whole foetal bovine serum. Upper and receiver 500ml. In polystyrene housing with polystyrene receiver connected by 1½ turn screw closure. Side-arm features a cellulosic vent plug and quick-disconnect tubing adapter. Colour coded yellow. Supplied sterile in cases of 12 units.

Cat. No.	Code	Pore size, µm
FU/36408	162-0020	0.2
FU/36418	162-0045	0.45

Filters, bottle top, SFCA

NALGENE

Tissue culture filters with 50, 75 or 90mm membranes. Securely screw onto media bottles with 33 or 45mm neck size. Allow direct filtration into sterile glass media bottles*. Moulded-in graduations. Side-arm features a quick-disconnect tubing adapter. Colour coded yellow. Supplied sterile in cases of 12 units.

Cat. No.	Code	Capacity ml	Pore size µm	Membrane dia. mm
For 33mm neck bottles				
FB/36302	290-3320	150	0.2	50
FB/36312	290-3345	150	0.45	50
FB/36304	291-3320	500	0.2	75
FB/36314	291-3345	500	0.45	75
FB/36316	292-3320	1000	0.2	90
For 45mm neck bottles				
FB/36402	290-4520	150	0.2	50
FB/36412	290-4545	150	0.45	50
FB/36404	291-4520	500	0.2	75
FB/36414	291-4545	500	0.45	75
FB/36416	292-4520	1000	0.2	90

*Caution: Use only sterile bottle designed for use with vacuum. Always use safety shield during vacuum procedures.

Filter units, PES membrane, screw closure

NALGENE

Ideal for biological and pharmaceutical sterilisation requirements. Graduations on both uppers and receivers. The side-arm has a cellulosic vent plug and quick-disconnect tubing adapter. In polystyrene housing connected to polystyrene receiver by leakproof threaded closure. Colour coded blue. Supplied sterile in cases of 12.

Cat. No.	Code	Upper/Receiver capacity	Membrane diameter
Pore size 0.45 µm			
FM/36212	165-0045	150 ml	50mm
FM/36214	168-0045	250 ml	50mm
FM/36216	166-0045	500 ml	75mm
FM/36217	169-0045	500 ml	90mm
FM/36218	167-0045	1000 ml	90mm